



# danumed® Extension Set for Gastrostomy Button

**Security  
meets Comfort**

**Durable TRITAN™  
Connector Material ✓**

**Kink-resistant  
Polyurethane Tube ✓**

**Innovative Rotating  
Ring ENSwivel® ✓**



**German  
Manufacturer ✓**

**Reliable  
Supply Chain ✓**

**High-End Button Connectors,  
Straight and Angled ✓**

**Long Duration of Use  
up to 14 Days ✓**


**Optimised  
Product Design**



# danumed® Extension Set

The danumed® Extension Set is designed for the administration of enteral nutrition, liquid medication, and decompression in connection with a Gastrostomy Button.

Innovative features ensure good flow rate and simplify the handling for the user.

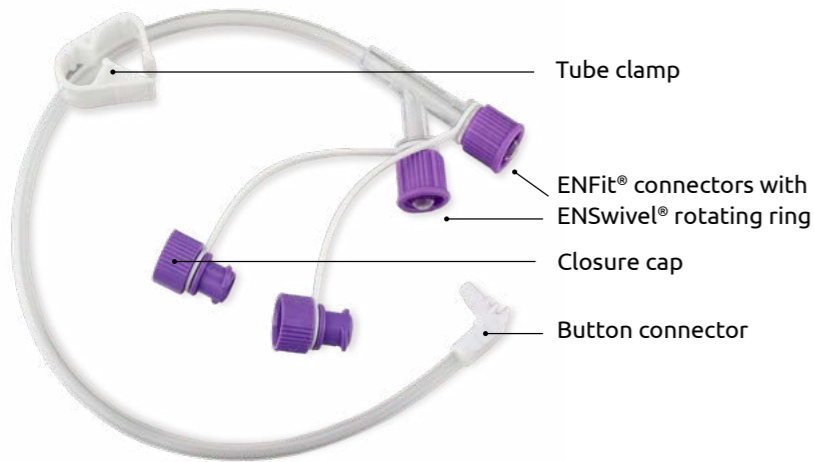

High-quality materials allow for up to 14 days of use:

### TRITAN™ Connector Material

- extremely resistant to aggressive drugs and solutions
- break-resistant, can avoid stress cracks

### Polyurethane Tube Material

- kink-resistant
- biocompatible, without harmful substances
- stays soft (WITHOUT plasticizers)


## Innovative Rotating Connector ENSwivel®

- For enhanced patient comfort and easy user handling
- Secure, tight connection with just a "Two-Finger-Turn". Simple, fast, comfortable.
- Twisting of tubes is avoided, as only the rotating ENSwivel® ring is turned and NOT the entire system.
- Fitted on all ENFit® connectors



## High-End Button Connector

- Rounded inner profile for unimpeded flow
- Easy, secure and tight connection to the Button
- Fits on all common gastrostomy buttons

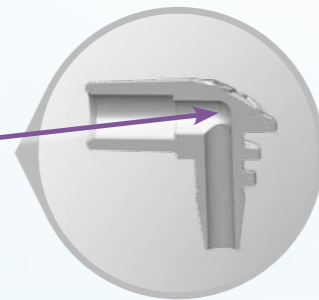

## Numerous Versions and Lengths


- Available ENFit® connector versions: single ENFit®, Y-ENFit®, ENFit® Bolus
- Available in different lengths: 5 cm, 30 cm, 60 cm



## Product Versions


### Bolus Nutrition with Syringe

ENFit® Bolus

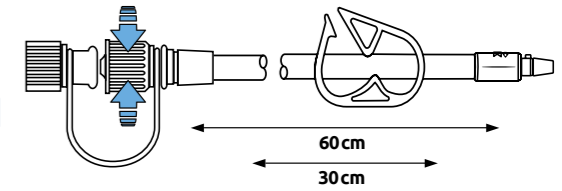


### Continuous Nutrition with Pump or Gravity Set

ENFit®

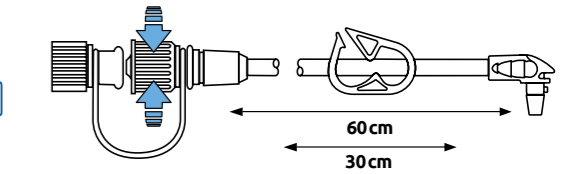

Y-ENFit®

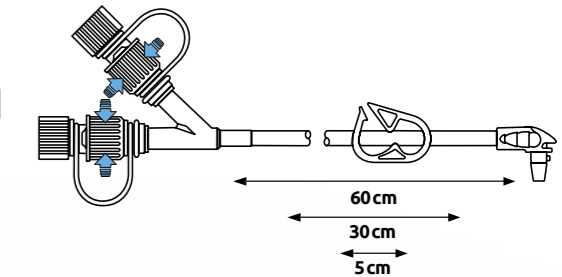

### danumed® Extension Set with angled Button Connector

REF	Description	Length	Feeding connector
DA 1010 3011	Extension Set <b>non-sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®
DA 1010 6011	Extension Set <b>non-sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®
DA 1010 3021	Extension Set <b>non-sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®
DA 1010 6021	Extension Set <b>non-sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®
DA 1011 3011	Extension Set <b>sterile</b> , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®
DA 1011 6011	Extension Set <b>sterile</b> , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®
DA 1011 0521	Extension Set <b>sterile</b> , 5 cm, Y-ENFit®, ENSwivel®	5 cm	Y-ENFit®
DA 1011 3021	Extension Set <b>sterile</b> , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®
DA 1011 6021	Extension Set <b>sterile</b> , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®

Delivered with right-angled button connector and **standard bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box.



### danumed® Extension Set with straight Button Connector

REF	Description	Length	Feeding connector
DA 1010 3032	Extension Set <b>non-sterile</b> , 30 cm, ENFit® Bolus, ENSwivel®	30 cm	ENFit® Bolus
DA 1010 6032	Extension Set <b>non-sterile</b> , 60 cm, ENFit® Bolus, ENSwivel®	60 cm	ENFit® Bolus
DA 1011 3032	Extension Set <b>sterile</b> , 30 cm, ENFit® Bolus, ENSwivel®	30 cm	ENFit® Bolus
DA 1011 6032	Extension Set <b>sterile</b> , 60 cm, ENFit® Bolus, ENSwivel®	60 cm	ENFit® Bolus

Delivered with straight button connector and **wide bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box. The ENFit® Bolus Extension Sets are available on request with funnel connector instead of ENFit® connector.



## **danumed Medizintechnik GmbH**

*Individual solutions are our standard.*

We are a family-run, innovative, German medical technology company with high quality standards, many years of expertise in enteral nutrition and a clear customer focus.

Individual customer requests are met with flexibility and reliability.  
Personal, competent, collaborative.

danumed Medizintechnik GmbH • Stadthof 18 • 93059 Regensburg • Germany  
Tel + 49 (0) 941/ 890 588-80 • Fax + 49 (0) 941/ 890 588-81

[info@danumed.com](mailto:info@danumed.com) • [www.danumed.com](http://www.danumed.com)