

INNOVATORS IN ENTERAL CARE



danumed
MEDIZINTECHNIK



Global Trust in German Design & Quality



danumed Medizintechnik GmbH

Individual solutions are our standard.

We are an owner-managed innovative German medical technology company with the highest quality standards, many years of experience in enteral nutrition and a clear customer orientation.

Our focus is on development, production and distribution of high-quality medical products for enteral nutrition. danumed offers a comprehensive product portfolio in enteral application technology for the safe and comfortable administration of enteral feeds.

All our customers appreciate our personal approach, our reliability and the fact that we always have an open ear for special requests.

Best regards from Regensburg, Germany,

Sonja and Hans Winklmann

Special product features of our danumed enteral application products:



ENFit®

Patient safety is a top priority for danumed. It is a matter of course that all our enteral products comply with the standards DIN EN ISO 20695 and 80369-3 and are equipped with ENFit® connectors, to avoid misconnections with other systems and to ensure the protection of enteral-fed patients worldwide.



ENSwivel®



danumed permanently improves its enteral products to enhance patient safety and comfort. A recent improvement is the equipment of most of our ENFit® connectors with the innovative rotating ring ENSwivel®, which makes the life of patients and users so much easier.

The ENSwivel® allows for easy connection and disconnection of ENFit® devices, even very tight connections can be released effortlessly. Furthermore, this feature avoids any twisting of the enteral tube, which is both painful for the patient and inconvenient for the user.

The danumed ENSwivel® – an innovative change with a big impact on patient comfort and user handling.



TRITAN™

danumed ENFit® connectors are made of TRITAN™, a light, robust and durable plastic material that is resistant to aggressive solutions and medicines. TRITAN™ ensures an extremely break resistant application and can largely avoid stress cracks.

All our enteral application products are manufactured with high-quality materials and are all latex-free, DEHP-free and BPA-free. The majority of our products is also PVC-free. Silicone and polyurethane are our preferred tube materials.

TABLE OF CONTENTS

danuButton®	
Compact XS / Compact Set / Premium Set.....	4
danumed® Extension Set for Gastrostomy Button.....	6
danumed® GastroTube	8
danumed® Enteral Syringes, Single Use	10
danumed® Enteral Syringes, Single-Patient Use	12
danumed® FeedBag.....	14
danumed® Gravity FeedBag.....	15
danumed® Gravity Spike and Gravity Universal.....	16
danumed® StomaMeter	18
danumed® Stopcock.....	19



danuButton® Compact XS / Compact Set / Premium Set

The danuButton® is a replacement feeding tube for long-term enteral nutrition and use with completely healed stoma. The low-profile feeding tube is easy to conceal under clothing and especially for active, mobile patients. The extensive accessories of the danuButton® Sets enable you to an easy and comfortable replacement of the gastrostomy probe, additional components are no longer necessary.

The decision for the danuButton® is also the decision for more quality of life through an innovative, high-quality product.



COMFORT

- **Large selection of sizes and lengths** for the individual and perfect fit
- **Extensive range of accessories** in the Sets for easy and convenient replacement of the gastrostomy tube – additional components not necessary
- **High comfort** due to soft, flexible silicone
- **Low-profile design** for more mobility and quality of life
- **Easy cleaning** and maintenance
- **Tilted** balloon valve

SAFETY


- Made of **high-quality medical grade silicone** with high clarity
- All sets are **latex-free, PVC-free, DEHP-free, BPA-free**
- **Radiopaque stripe** along the entire shaft
- **Apple-shaped inner retention balloon with high balloon stability** for long usage, secure fixation and to avoid leakage
- **Flat, rounded tip** enables atraumatic insertion
- **Recessed tip** to avoid trauma
- Secure **anti-reflux** valve

The danumed® Extension Sets included in the danuButton® Compact and Premium Sets are made of kink-resistant, biocompatible polyurethane. They are equipped with ENFit® connectors according to DIN EN ISO 80369-3 and danumed's innovative rotating connector, the ENSwivel®.



The danuButton® Compact Set contains:

- danuButton® Gastrostomy Button
- danumed® Extension Set for Gastrostomy Button, 30 cm, Y-ENFit®, ENSwivel®, angled button connector
- danumed® Introducer
- danumed® Balloon Syringe 10 ml
- optiLube lubricant
- Gauze Pads (4 x)

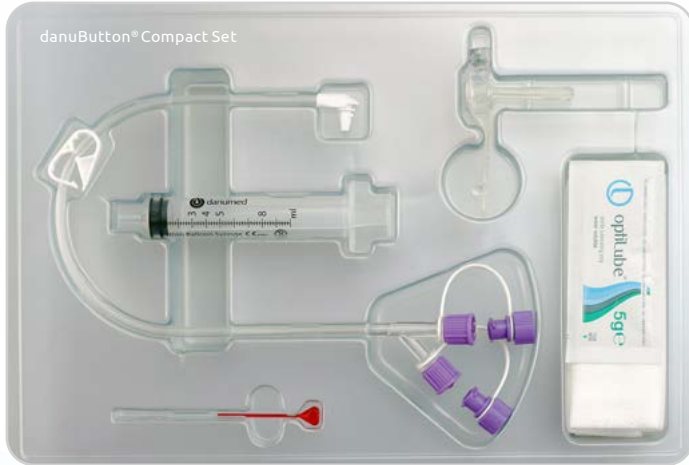

		Fr-Size					
		12	14	16	18	20	24
Length (cm)	0.8	DA 1042 1208	DA 1042 1408	DA 1042 1608	DA 1042 1808	DA 1042 2008	DA 1042 2408
	1.0	DA 1042 1210	DA 1042 1410	DA 1042 1610	DA 1042 1810	DA 1042 2010	DA 1042 2410
	1.2	DA 1042 1212	DA 1042 1412	DA 1042 1612	DA 1042 1812	DA 1042 2012	DA 1042 2412
	1.5	DA 1042 1215	DA 1042 1415	DA 1042 1615	DA 1042 1815	DA 1042 2015	DA 1042 2415
	1.7	DA 1042 1217	DA 1042 1417	DA 1042 1617	DA 1042 1817	DA 1042 2017	DA 1042 2417
	2.0	DA 1042 1220	DA 1042 1420	DA 1042 1620	DA 1042 1820	DA 1042 2020	DA 1042 2420
	2.3	DA 1042 1223	DA 1042 1423	DA 1042 1623	DA 1042 1823	DA 1042 2023	DA 1042 2423
	2.5	DA 1042 1225	DA 1042 1425	DA 1042 1625	DA 1042 1825	DA 1042 2025	DA 1042 2425
	2.7	DA 1042 1227	DA 1042 1427	DA 1042 1627	DA 1042 1827	DA 1042 2027	DA 1042 2427
	3.0	DA 1042 1230	DA 1042 1430	DA 1042 1630	DA 1042 1830	DA 1042 2030	DA 1042 2430
	3.5	DA 1042 1235	DA 1042 1435	DA 1042 1635	DA 1042 1835	DA 1042 2035	DA 1042 2435
	4.0	DA 1042 1240	DA 1042 1440	DA 1042 1640	DA 1042 1840	DA 1042 2040	DA 1042 2440
	4.5	DA 1042 1245	DA 1042 1445	DA 1042 1645	DA 1042 1845	DA 1042 2045	DA 1042 2445
	5.0	DA 1042 1250	DA 1042 1450	DA 1042 1650	DA 1042 1850	DA 1042 2050	DA 1042 2450
	5.5		DA 1042 1455				
	6.0				DA 1042 1860		
6.5		DA 1042 1465					

1 Set per box, 5 boxes per shipping box.

The danuButton® Compact XS contains:


- danuButton® Gastrostomy Button
- danumed® Introducer
- danumed® Balloon Syringe 10 ml

		Fr-Size					
		12	14	16	18	20	24
Length (cm)	0.8	DA 1043 1208	DA 1043 1408	DA 1043 1608	DA 1043 1808	DA 1043 2008	DA 1043 2408
	1.0	DA 1043 1210	DA 1043 1410	DA 1043 1610	DA 1043 1810	DA 1043 2010	DA 1043 2410
	1.2	DA 1043 1212	DA 1043 1412	DA 1043 1612	DA 1043 1812	DA 1043 2012	DA 1043 2412
	1.5	DA 1043 1215	DA 1043 1415	DA 1043 1615	DA 1043 1815	DA 1043 2015	DA 1043 2415
	1.7	DA 1043 1217	DA 1043 1417	DA 1043 1617	DA 1043 1817	DA 1043 2017	DA 1043 2417
	2.0	DA 1043 1220	DA 1043 1420	DA 1043 1620	DA 1043 1820	DA 1043 2020	DA 1043 2420
	2.3	DA 1043 1223	DA 1043 1423	DA 1043 1623	DA 1043 1823	DA 1043 2023	DA 1043 2423
	2.5	DA 1043 1225	DA 1043 1425	DA 1043 1625	DA 1043 1825	DA 1043 2025	DA 1043 2425
	2.7	DA 1043 1227	DA 1043 1427	DA 1043 1627	DA 1043 1827	DA 1043 2027	DA 1043 2427
	3.0	DA 1043 1230	DA 1043 1430	DA 1043 1630	DA 1043 1830	DA 1043 2030	DA 1043 2430
	3.5	DA 1043 1235	DA 1043 1435	DA 1043 1635	DA 1043 1835	DA 1043 2035	DA 1043 2435
	4.0	DA 1043 1240	DA 1043 1440	DA 1043 1640	DA 1043 1840	DA 1043 2040	DA 1043 2440
	4.5	DA 1043 1245	DA 1043 1445	DA 1043 1645	DA 1043 1845	DA 1043 2045	DA 1043 2445
	5.0	DA 1043 1250	DA 1043 1450	DA 1043 1650	DA 1043 1850	DA 1043 2050	DA 1043 2450
	5.5		DA 1043 1455				
	6.0		DA 1043 1460		DA 1043 1860		
6.5		DA 1043 1465					

1 Set per box, 5 boxes per shipping box.

The danuButton® Premium Set contains:

- danuButton® Gastrostomy Button
- danumed® Extension Set for Gastrostomy Button, 30 cm, Bolus-ENFit®, ENSwivel®, straight button connector
- danumed® Extension Set for Gastrostomy Button, 30 cm, Y-ENFit®, ENSwivel®, angled button connector
- danumed® Introducer
- danumed® Balloon Syringe 10 ml
- danumed® Enteral Syringe 60 ml, ENFit®
- optiLube lubricant
- Gauze Pads (4 x)

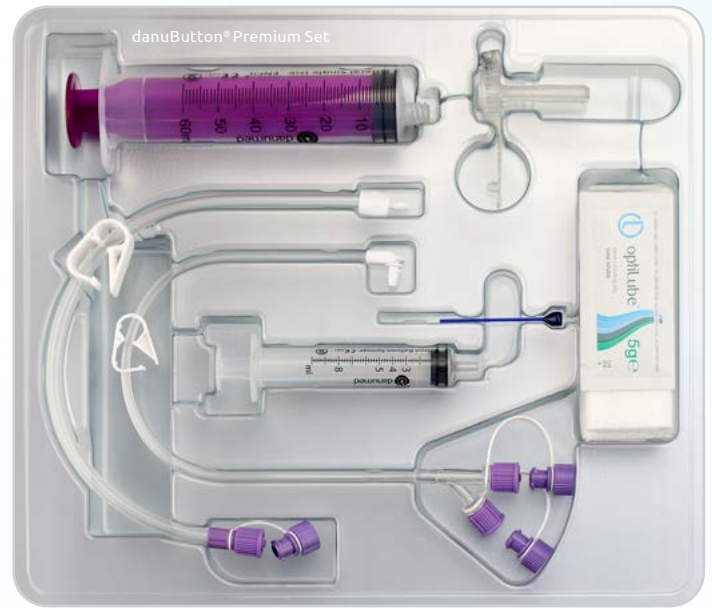

		Fr-Size					
		12	14	16	18	20	24
Length (cm)	0.8	DA 1041 1208	DA 1041 1408	DA 1041 1608	DA 1041 1808	DA 1041 2008	DA 1041 2408
	1.0	DA 1041 1210	DA 1041 1410	DA 1041 1610	DA 1041 1810	DA 1041 2010	DA 1041 2410
	1.2	DA 1041 1212	DA 1041 1412	DA 1041 1612	DA 1041 1812	DA 1041 2012	DA 1041 2412
	1.5	DA 1041 1215	DA 1041 1415	DA 1041 1615	DA 1041 1815	DA 1041 2015	DA 1041 2415
	1.7	DA 1041 1217	DA 1041 1417	DA 1041 1617	DA 1041 1817	DA 1041 2017	DA 1041 2417
	2.0	DA 1041 1220	DA 1041 1420	DA 1041 1620	DA 1041 1820	DA 1041 2020	DA 1041 2420
	2.3	DA 1041 1223	DA 1041 1423	DA 1041 1623	DA 1041 1823	DA 1041 2023	DA 1041 2423
	2.5	DA 1041 1225	DA 1041 1425	DA 1041 1625	DA 1041 1825	DA 1041 2025	DA 1041 2425
	2.7	DA 1041 1227	DA 1041 1427	DA 1041 1627	DA 1041 1827	DA 1041 2027	DA 1041 2427
	3.0	DA 1041 1230	DA 1041 1430	DA 1041 1630	DA 1041 1830	DA 1041 2030	DA 1041 2430
	3.5	DA 1041 1235	DA 1041 1435	DA 1041 1635	DA 1041 1835	DA 1041 2035	DA 1041 2435
	4.0	DA 1041 1240	DA 1041 1440	DA 1041 1640	DA 1041 1840	DA 1041 2040	DA 1041 2440
	4.5	DA 1041 1245	DA 1041 1445	DA 1041 1645	DA 1041 1845	DA 1041 2045	DA 1041 2445
	5.0	DA 1041 1250	DA 1041 1450	DA 1041 1650	DA 1041 1850	DA 1041 2050	DA 1041 2450
	5.5		DA 1041 1455				
	6.0		DA 1041 1460		DA 1041 1860		
6.5		DA 1041 1465					

1 Set per box, 5 boxes per shipping box.



All danuButton® Sets are:

- Latex-free
- PVC-free
- DEHP-free
- BPA-free



danumed® Extension Set for Gastrostomy Button

Extension set for delivery of enteral nutrition, liquid medication and decompression.

Innovative features ensure good flow and easy handling.

Cost-efficient solution with usage for up to 14 days.

Optimised Product Design



COMFORT



Innovative Rotating Connector ENSwivel

- Allows for the easy connection and release of ENFit® connections with a **"Two-Finger-Turn"**. Simple, fast and comfortable.
- Only the rotating Swivel is turned and not the entire connector. This feature **avoids any twisting of the enteral tube** which can be painful for the patient and inconvenient for the user.
- All danumed ENFit® connectors offer the advantages of the ENSwivel®: **ENFit®, Y-ENFit®, ENFit® Bolus**.

High-End Button Connector

- The design of the locking mechanism ensures **easy, safe and tight connection** to the button.



Rounded kink in the inner lumen for unimpeded flow.



Ergonomic form of the button connector with ribbed profile for good grip.

SAFETY

Maximum safety for the patient

- The ENFit® connector according to DIN EN ISO 80369-3 **avoids tubing misconnections** and enhances patient safety significantly.
- All ENFit® connectors from danumed are made of **TRITAN™**, a plastic material resistant to aggressive solutions/medicines.
- All extension sets with ENFit® connector are **free from latex, PVC, DEHP and BPA**.
- The tube is made of polyurethane – a biocompatible and kink-resistant material that allows duration of use for **up to 14 days**.



danumed® Extension Set with angled Button Connector



REF	Description	Length	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-free
DA 1010 3011	Extension Set non-sterile , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®	●	●	●	●
DA 1010 6011	Extension Set non-sterile , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®	●	●	●	●
DA 1010 3021	Extension Set non-sterile , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®	●	●	●	●
DA 1010 6021	Extension Set non-sterile , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®	●	●	●	●
DA 1011 3011	Extension Set, sterile , 30 cm, ENFit®, ENSwivel®	30 cm	ENFit®	●	●	●	●
DA 1011 6011	Extension Set sterile , 60 cm, ENFit®, ENSwivel®	60 cm	ENFit®	●	●	●	●
DA 1011 0521	Extension Set sterile , 5 cm, Y-ENFit®, ENSwivel®	5 cm	Y-ENFit®	●	●	●	●
DA 1011 3021	Extension Set sterile , 30 cm, Y-ENFit®, ENSwivel®	30 cm	Y-ENFit®	●	●	●	●
DA 1011 6021	Extension Set sterile , 60 cm, Y-ENFit®, ENSwivel®	60 cm	Y-ENFit®	●	●	●	●

Delivered with right-angled button connector and **standard bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box.



danumed® Extension Set with straight Button Connector



REF	Description	Length	Feeding connector	Latex-free	PVC-free	DEHP-free	BPA-free
DA 1010 3032	Extension Set non-sterile , 30 cm, ENFit® Bolus, ENSwivel®	30 cm	ENFit® Bolus	●	●	●	●
DA 1010 6032	Extension Set non-sterile , 60 cm, ENFit® Bolus, ENSwivel®	60 cm	ENFit® Bolus	●	●	●	●
DA 1011 3032	Extension Set sterile , 30 cm, ENFit® Bolus®, ENSwivel®	30 cm	ENFit® Bolus	●	●	●	●
DA 1011 6032	Extension Set sterile , 60 cm, ENFit® Bolus, ENSwivel®	60 cm	ENFit® Bolus	●	●	●	●

Delivered with straight button connector and **wide bore tube**. Individually packaged. 10 pieces in a box, 10 boxes in a shipping box. The Bolus Extension Sets are available on request with funnel connector instead of ENFit® connector.

Related Products to danumed® Extension Sets:



danuButton®
Gastrostomy Button



danumed®
Enteral Syringes, Single Use
Enteral Syringes, Single-Patient Use



danumed® GastroTube



Optimised Product Design



ENSwivel



90° danuPlate



danuClip*

Gastrostomy feeding tube with retention balloon for long-term enteral nutrition

For superior patient comfort and ease of use:

- **Innovative rotating ring ENSwivel®**, for easy connection and release of ENFit® connections, without twisting the tubes. Secure, tight connection with just a "Two-Finger-Turn". Simple, fast, comfortable.
- **90° danuPlate** – soft external retention plate with 90° angulation for secure and discreet fixation.
- **danuClip** – tube clamp, prevents reflux of nutrition/stomach contents*.

*Available as a separate product, please see REF on back page.

COMFORT

- **Excellent patient comfort** due to soft, flexible silicone
- **Rotating ring ENSwivel®** on all ENFit® connectors
- **90° danuPlate** – structured design allows for air circulation in the stoma site area
- **danuClip** – practical, simple and quick to use*
- **Extensive selection of sizes** from 12 Fr to 24 Fr for the individual and perfect fit
- Balloon fill valve with Luer connector, with **colour coding** according to Fr size
- **Length marking** in cm along the tube
- **Display of the balloon volume** on the shaft
- **Convenient accessories** with every GastroTube: a pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10 ml** to inflate the retention balloon



SAFETY

- Gastric replacement feeding tube made of **high-quality medical grade silicone**
- **Free of Latex, PVC, DEHP and BPA**
- **Flat rounded tip** to reduce trauma risk
- **Inner retention balloon with high balloon stability** for long usage, secure fixation and to avoid leakage
- **ENFit® connectors** according to DIN EN ISO 80369-3 to prevent tubing misconnections with other systems
- Available with **ENFit® and Y-ENFit®** connectors
- **Radiopaque stripe** along the shaft

Easy handling and secure, discreet fixation of the GastroTube with the 90° danuPlate:



danumed® GastroTubes with ENFit® and Y-ENFit® Supplied in a set with lubricant and danumed® balloon syringe 10 ml.



REF	Product description	Size in Fr	Pieces/sales unit	Latex-free	PVC-free	DEHP-free	Bisphenol-A-free
DA 1044 0112	danumed® GastroTube with ENFit® and ENSwivel®	12	1 piece	●	●	●	●
DA 1044 0114	danumed® GastroTube with ENFit® and ENSwivel®	14	1 piece	●	●	●	●
DA 1044 0116	danumed® GastroTube with ENFit® and ENSwivel®	16	1 piece	●	●	●	●
DA 1044 0118	danumed® GastroTube with ENFit® and ENSwivel®	18	1 piece	●	●	●	●
DA 1044 0120	danumed® GastroTube with ENFit® and ENSwivel®	20	1 piece	●	●	●	●
DA 1044 0122	danumed® GastroTube with ENFit® and ENSwivel®	22	1 piece	●	●	●	●
DA 1044 0124	danumed® GastroTube with ENFit® and ENSwivel®	24	1 piece	●	●	●	●
DA 1044 0212	danumed® GastroTube with Y-ENFit® and ENSwivel®	12	1 piece	●	●	●	●
DA 1044 0214	danumed® GastroTube with Y-ENFit® and ENSwivel®	14	1 piece	●	●	●	●
DA 1044 0216	danumed® GastroTube with Y-ENFit® and ENSwivel®	16	1 piece	●	●	●	●
DA 1044 0218	danumed® GastroTube with Y-ENFit® and ENSwivel®	18	1 piece	●	●	●	●
DA 1044 0220	danumed® GastroTube with Y-ENFit® and ENSwivel®	20	1 piece	●	●	●	●
DA 1044 0222	danumed® GastroTube with Y-ENFit® and ENSwivel®	22	1 piece	●	●	●	●
DA 1044 0224	danumed® GastroTube with Y-ENFit® and ENSwivel®	24	1 piece	●	●	●	●
DA 1044 0700	danuClip Tube Clamp	--	15 pcs.	●	●	●	●

Each danumed® GastroTube is individually sterile packed. The danumed® GastroTube is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10 ml in each single box. 10 boxes in a shipping carton.

Related products to danumed® GastroTubes:



danumed®
Enteral Syringes,
Single Use and
Single-Patient Use



danumed® Gravity Spike,
Gravity Universal and
Gravity Feedbag
Gravity feeding sets



danumed® Enteral Syringes, Single Use

danumed® Enteral Syringes, Single Use, with female ENFit® connector according to DIN EN ISO 80369-3

For the administration of enteral feeds, medicines and liquids, for decompression and flushing

Extensive range of sizes for safe use in neonates, children and adults



COMFORT

Comfortable handling for the user

- **Transparent syringe cylinder** for viewing the content
- **Precisely legible** and wipe-proof scaling
- **Secure hold** with ribbed grip plate
- **Smooth**, even plunger movement
- **Noticeable safety stop** to prevent unintentional removal of the plunger from the syringe cylinder

danumed® Syringe Cap

One size fits all:
Universal syringe cap for all sizes and versions of danumed® Enteral Syringes.



danumed® Straw

Round-tip enteral straw for safe and clean withdrawal of enteral nutrition, medication and liquids.



With ENFit® connector.

Available in different sizes and lengths.

SAFETY

Excellent safety for the patient

- **ENFit® connector** to reduce risk of tubing misconnections
- **"Low-Dose" Tip** in 1 ml, 3 ml and 5 ml syringes
- High quality biocompatible material **free of latex, PVC, DEHP and BPA**
- **Validated grading** for high measurement accuracy



danumed® Enteral Syringes, Single Use, are individually sterile packed and available in standard sizes of 1 ml, 3 ml, 5 ml, 10 ml, 20 ml, 60 ml and 100 ml. Further sizes are available on request.

Safe products according to the standards:

DIN EN ISO 80369-3, DIN EN ISO 7886-1, DIN EN ISO 20695

danumed® Enteral Syringes, Single Use

REF	Product description	Size	Pieces/ sales unit	Latex- free	PVC- free	DEHP- free	Bisphenol- A-free
DA 1021 0110	danumed® Enteral Syringe, single use, with centric tip, low-dose tip and ENFit®	1 ml	100	●	●	●	●
DA 1021 0310	danumed® Enteral Syringe, single use, with centric tip, low-dose tip and ENFit®	3 ml	100	●	●	●	●
DA 1021 0510	danumed® Enteral Syringe, single use, with centric tip, low-dose tip and ENFit®	5 ml	100	●	●	●	●
DA 1021 0500	danumed® Enteral Syringe, single use, with centric tip and ENFit®	5 ml	100	●	●	●	●
DA 1021 1000	danumed® Enteral Syringe, single use, with centric tip and ENFit®	10 ml	100	●	●	●	●
DA 1021 2000	danumed® Enteral Syringe, single use, with centric tip and ENFit®	20 ml	100	●	●	●	●
DA 1021 6001	danumed® Enteral Syringe, single use, with eccentric tip and ENFit®	60 ml	100	●	●	●	●
DA 1021 9900	danumed® Enteral Syringe, single use, with centric tip and ENFit®	100 ml	25	●	●	●	●

All products are disposable and individually sterile packed. Further sizes are available on request.

danumed® Syringe Caps and Enteral Straws

REF	Product description	Size/ length	Pieces/ sales unit	Latex- free	PVC- free	DEHP- free	Bisphenol- A-free
DA 1021 0055	danumed® Syringe Caps for Enteral Syringes, single use, ENFit®	one size	100	●	●	●	●
DA 1021 1210	danumed® Enteral Straw for milk, ENFit®	12 Fr/ 100 mm	30	●	●	●	●
DA 1021 1215	danumed® Enteral Straw for milk, ENFit®	12 Fr/ 150 mm	30	●	●	●	●
DA 1021 1604	danumed® Enteral Straw for milk, ENFit®	16 Fr/ 40 mm	30	●	●	●	●
DA 1021 1615	danumed® Enteral Straw for milk, ENFit®	16 Fr/ 150 mm	30	●	●	●	●
DA 1021 0605	danumed® Enteral Straw for medication, ENFit®	6 Fr/ 50 mm	30	●	●	●	●

All products are disposable and individually sterile packed.

Related products to danumed® Enteral Syringes, Single Use:



danumed®
Syringe Cap



danumed®
Straw



danumed® Extension Set
for Gastrostomy Button



danuButton®
Gastrostomy Button



danumed® Enteral Syringes, Single-Patient Use

danumed® Enteral Syringes, Single-Patient Use, with female ENFit® connector according to EN ISO 80369-3

For the administration of enteral feeds, medicines and liquids, for decompression and flushing

Extensive range of sizes for safe use in neonates, children and adults

Cost-efficient solution with usage for up to 14 days in home use



COMFORT

Comfortable handling for the user

- **Washable and reusable up to 14 days**
- **Transparent syringe cylinder** for viewing the content
- **Precisely legible** and wipe-proof scaling
- **Secure hold** with ribbed grip plate
- **Smooth**, even plunger movement
- **Noticeable safety stop** to prevent unintentional removal of the plunger from the syringe cylinder

danumed® Syringe Cap

One size fits all:
Universal syringe cap for all sizes and versions of danumed® Enteral Syringes.



danumed® Straw

Round-tip enteral straw for safe and clean withdrawal of enteral nutrition, medication and liquids.

With ENFIT® connector.

Available in different sizes and lengths.



SAFETY

Excellent safety for the patient

- **ENFit® connector** to reduce risk of tubing misconnections
- **“Low-Dose” Tip** in 1 ml and 3 ml syringes
- High quality biocompatible material **free of latex, PVC, DEHP and BPA**
- **Validated grading** for high measurement accuracy



danumed® Enteral Syringes, Single Patient Use, are individually sterile packed and available in standard sizes of 1 ml, 3 ml, 5 ml, 10 ml, 20 ml and 60 ml. Further sizes are available on request.

Safe products according to the standards:
DIN EN ISO 80369-3, DIN EN ISO 7886-1,
DIN EN ISO 20695

danumed® Enteral Syringes, Single-Patient Use

REF	Product description	Size	Pieces/ sales unit	Latex- free	PVC-free	DEHP-free	Bisphenol- A-free
DA 1022 0110	danumed® Enteral Syringe, single-patient use, with centric tip, low-dose tip and ENFit®	1 ml	100	●	●	●	●
DA 1022 0310	danumed® Enteral Syringe, single-patient use, with centric tip, low-dose tip and ENFit®	3 ml	100	●	●	●	●
DA 1022 0500	danumed® Enteral Syringe, single-patient use, with centric tip and ENFit®	5 ml	100	●	●	●	●
DA 1022 1000	danumed® Enteral Syringe, single-patient use, with centric tip and ENFit®	10 ml	100	●	●	●	●
DA 1022 2000	danumed® Enteral Syringe, single-patient use, with centric tip and ENFit®	20 ml	100	●	●	●	●
DA 1022 6001	danumed® Enteral Syringe, single-patient use, with eccentric tip and ENFit®	60 ml	50	●	●	●	●

All products are reusable and individually sterile packed. Further sizes are available on request.

danumed® Syringe Caps and Enteral Straws

REF	Product description	Size/ length	Pieces/ sales unit	Latex- free	PVC-free	DEHP-free	Bisphenol- A-free
DA 1021 0055	danumed® Syringe Caps for Enteral Syringes, ENFit®	one size	100	●	●	●	●
DA 1021 1210	danumed® Enteral Straw for milk, ENFit®	12 Fr/ 100 mm	30	●	●	●	●
DA 1021 1215	danumed® Enteral Straw for milk, ENFit®	12 Fr/ 150 mm	30	●	●	●	●
DA 1021 1604	danumed® Enteral Straw for milk, ENFit®	16 Fr/ 40 mm	30	●	●	●	●
DA 1021 1615	danumed® Enteral Straw for milk, ENFit®	16 Fr/ 150 mm	30	●	●	●	●
DA 1021 0605	danumed® Enteral Straw for medication, ENFit®	6 Fr/ 50 mm	30	●	●	●	●

All products are disposable and individually sterile packed.

Related products to danumed® Enteral Syringes, Single-Patient Use:



danumed®
Syringe Cap



danumed®
Straw



danumed® Extension Set
for Gastrostomy Button



danuButton®
Gastrostomy Button



Global Trust in German Design & Quality



danumed® FeedBag

Empty bag for enteral nutrition and fluid balancing

The danumed® Feedbag is an empty feeding bag with a large filling volume of 1500 ml.

The ENPlus port on the outlet tube allows connection to the danumed® Gravity Sets and to all other commercially available enteral gravity sets with ENPlus connection.



SAFETY & COMFORT

- **ENPlus port** for connection to enteral feeding systems
- **Large bag volume** of 1500 ml
- **Inscribable bag material** for clear patient identification
- **Scaling** in 100ml-steps
- With tube clamp for **quick application stop**
- Large topfiller for **convenient filling** of the bag

danumed® FeedBag

The high-quality, economical version made of latex-free, DEHP- and BPA-free bag material.

REF	Product description	Pieces/sales unit
DA 1051 1000	danumed®FeedBag	30 pieces

The danumed®FeedBag is individually sterile packed. 30 pieces in a shipping box.

The danumed®FeedBags are

- Latex-free • DEHP-free • BPA-free

danumed® FeedBag Premium

The first-class PVC-free version made of EVA bag material and with premium topfiller.

REF	Product description	Pieces/sales unit
DA 1051 1001	danumed®FeedBag Premium	30 pieces

The danumed®FeedBag Premium is individually sterile packed. 30 pieces in a shipping box.

The danumed®FeedBags Premium are

- Latex-free • PVC-free • DEHP-free • BPA-free



danumed® Gravity FeedBag

Gravity Feeding Set – Enteral Bag 1500 ml with firmly connected tube system

For the administration of enteral nutrition, liquids and medication

COMFORT

- **Large bag volume** of 1500 ml
- With large topfiller for **convenient filling** of the bag
- **Large tube lumen** for good flow rate
- **Three-way stopcock** for easy and comfortable administration of medication
- Transparent **drip chamber**
- **Roller clamp** for flow rate adjustment
- **Inscribable bag material** for clear patient identification
- **Scaling** in 100ml-steps

SAFETY

- **High quality material** free of latex, DEHP and BPA
- **ENFit®** connectors in accordance with DIN EN ISO 80369-3 to reduce the risk of tubing misconnections and enhance patient safety
- **Protective caps** at the ENFit® connectors



REF	Product description	Pieces/sales unit
DA 1051 0210	danumed® GravityFeedBag	30 pieces

Each danumed® Gravity FeedBag is individually sterile packed. 30 pieces in a shipping box.

danumed® Gravity FeedBags are

- Latex-free
- DEHP-free
- BPA-free





danumed® Gravity Spike and Gravity Universal

Gravity Feeding Sets for the administration of enteral nutrition

Two versions for flexible connection possibilities to empty bags, prefilled feeding bags, crown caps and wide-neck bottles:

- danumed® Gravity Spike with ENPlus Spike
- danumed® Gravity Universal with ENPlus Spike and additional universal adapter

The duration of use of the danumed® Gravity Feeding Sets is up to 24 hours.



COMFORT

- **Flexibility** through the possibility of connecting to different food containers
- **Large tube lumen** for good flow rate
- **Three-way stopcock** for easy and comfortable administration of medication
- Transparent **drip chamber**
- **Roller clamp** for flow rate adjustment

SAFETY

- **ENFit® connectors** in accordance with DIN EN ISO 80369-3, to reduce the risk of tubing misconnections with other systems
- **Protective caps** on the ENPlus Spike and on all ENFit® connectors
- All danumed® Gravity Feeding Sets are **free of Latex, DEHP and BPA**

danumed® Gravity Spike

With **ENPlus Spike** for connection to pre-filled feeding bags and to the danumed®FeedBags.

REF	Product description	Pieces/sales unit
DA 1051 0010	danumed®Gravity Spike	30 pieces

The danumed®Gravity Spike is individually sterile packed. 30 pieces in a shipping box.



All danumed®Gravity Feeding Sets are
 • Latex-free • DEHP-free • BPA-free

danumed® Gravity Universal

With **ENPlus Spike** for connection to danumed®FeedBags or pre-filled feeding bags and with additional **Universal Adapter** with vent valve for connection to crown caps and wide-neck bottles.

REF	Product description	Pieces/sales unit
DA 1051 0110	danumed®Gravity Universal	30 pieces

The danumed®Gravity Universal is individually sterile packed. 30 pieces in a shipping box.



The danumed®Universal Adapter with vent valve is also available as a separate product:

REF	Product description	Pieces/sales unit
DA 1051 8003	danumed® Universal Adapter	30 pieces

Related products to danumed® Gravity Spike/Universal:



danumed®FeedBag



danumed®Universal Adapter



danumed®GastroTube






danumed® StomaMeter

Stoma Measuring Device with retention balloon

The StomaMeter is a measuring instrument for determining the stoma length. The length of the stoma tract is measured simply, accurately and reliably to ensure the perfect selection of the optimal gastrostomy button.

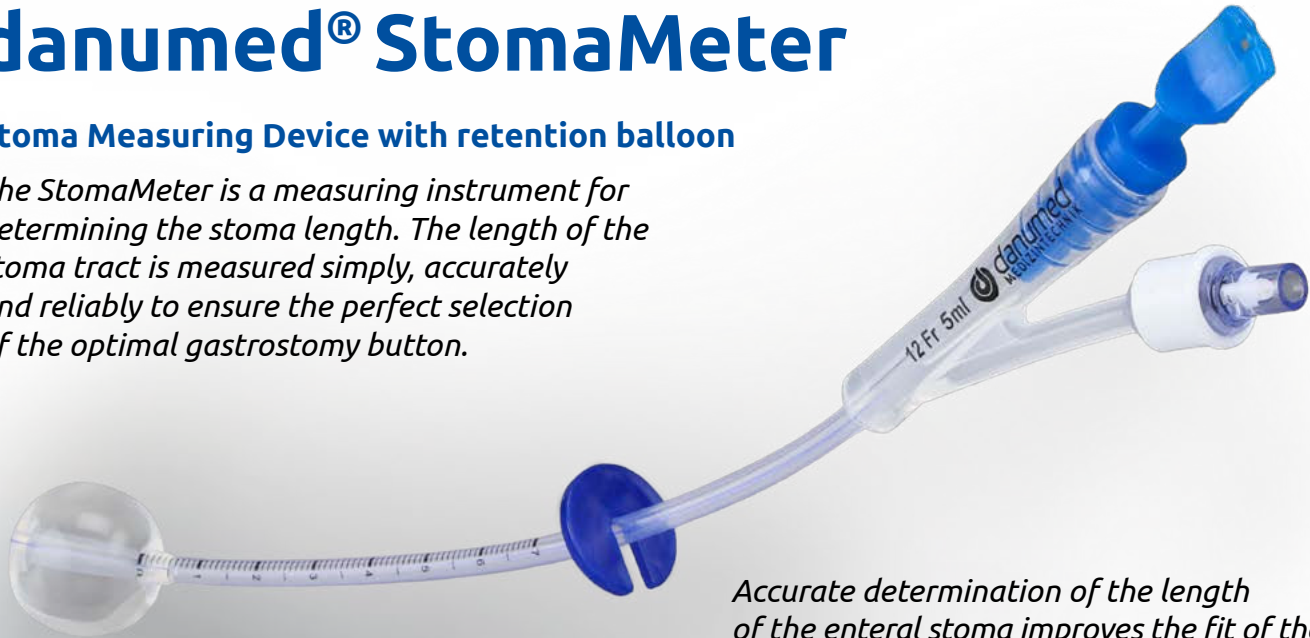

Accurate determination of the length of the enteral stoma improves the fit of the gastrostomy button and thus patient safety and comfort.

COMFORT

- **Soft, flexible material** provides patient comfort
- **Flat retention plate** simplifies reading of measured values
- **Imprint of the balloon volume** on the funnel
- **Convenient accessories** in every single package: A pack of **lubricant** for easy insertion and a **danumed® balloon syringe 10ml** to inflate the retention balloon.
- Available in **size 12 Fr**, suitable for each stoma



SAFETY

- Made of high-quality **medical grade silicone**
- **Latex-free, PVC-free, DEHP-free** and **BPA-free** material
- **Flat rounded tip** to avoid trauma
- **Radiopaque stripe** along the shaft
- **Precisely legible** and wipe-resistant scaling

danumed® StomaMeter

REF	Size in French	Pieces/sales unit
DA 1045 0112	12	1 piece

The danumed® StomaMeter stoma measuring device is individually sterile packed. The StomaMeter is delivered in a sterile peel-pouch and additionally with lubricant and a sterile danumed® balloon syringe 10ml in each single box. 10 boxes in a shipping carton.

The danumed® StomaMeter is

- **Latex-free** • **PVC-free** • **DEHP-free** • **BPA-free**



Global Trust in German Design & Quality



danumed® Stopcock

Three-way stopcock with ENFit® and ENSwivel® for enteral application systems

Patient safety and user comfort are top priorities for danumed. All of our enteral nutrition systems are equipped with ENFit® connectors in accordance with DIN EN ISO 80369-3 in order to reduce the risk of unintentional faulty connections with other systems.

danumed® Stopcock is an ENFit® three-way stopcock for the separate administration of enteral nutrition and medication. The male ENFit® connectors are equipped with the innovative rotating ring ENSwivel® for better patient comfort and easy handling.

Lightweight design for comfortable use.



COMFORT

ENSwivel®

- Twisting of the enteral tubes is avoided since not the whole connector is turned, but only the ENSwivel®
- Secure, tight connection with just a "Two-Finger-Turn". Simple, fast, comfortable.



Larger flow hole

- For consistent flow rate of enteral feeds, liquids and medication
- Opening in the form of a rounded bottom slot for unimpeded flow



SAFETY

ENFit®

- A female ENFit® connector for connection to a feeding tube
- Two male ENFit® connectors with ENSwivel® to expand connection options

TRITAN™

- The base unit of the danumed® Stopcock is made of TRITAN™ – a light, robust and durable plastic that is resistant to aggressive solutions and medication
- Stress cracks are largely avoided, the danumed® Stopcock is extremely break-resistant in the application

danumed® Stopcock

REF	Product description	Piece/sales unit
DA 1051 7022	danumed® Stopcock	20 pieces

Each danumed® Stopcock is separately sterile-packed. 20 pieces in each box, 800 pieces in each shipping box.

The danumed® Stopcock is

- Latex-free • PVC-free • DEHP-free • BPA-free



danumed Medizintechnik GmbH

Individual solutions are our standard.

We are an innovative, German medical technology company with the highest quality standards, many years of expertise in enteral nutrition and clear customer focus. We also meet individual customer requests with flexibility and reliability. Personal, competent, collaborative.

danumed Medizintechnik GmbH • Stadthof 18 • 93059 Regensburg • Germany
Tel + 49 (0) 941/ 890 588-80 • Fax + 49 (0) 941/ 890 588-81
info@danumed.com • www.danumed.com