




# danumed® Gravity Spike and Gravity Universal

## Gravity Feeding Sets for the administration of enteral nutrition

Two versions for flexible connection possibilities to empty bags, prefilled feeding bags, crown caps and wide-neck bottles:

- danumed® Gravity Spike with ENPlus Spike
- danumed® Gravity Universal with ENPlus Spike and additional universal adapter

The duration of use of the danumed® Gravity Feeding Sets is up to 24 hours.

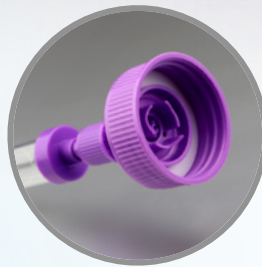

### COMFORT

- **Flexibility** through the possibility of connecting to different food containers
- **Large tube lumen** for good flow rate
- **Three-way stopcock** for easy and comfortable administration of medication
- Transparent **drip chamber**
- **Roller clamp** for flow rate adjustment

### SAFETY

- **ENFit® connectors** in accordance with DIN EN ISO 80369-3, to reduce the risk of tubing misconnections with other systems
- **Protective caps** on the ENPlus Spike and on all ENFit® connectors
- All danumed® Gravity Feeding Sets are **free of Latex, DEHP and BPA**

## danumed® Gravity Spike

With **ENPlus Spike** for connection to pre-filled feeding bags and to the danumed® FeedBags.

REF	Product description	Pieces/sales unit
DA 1051 0010	danumed® Gravity Spike	30 pieces

The danumed® Gravity Spike is individually sterile packed. 30 pieces in a shipping box.



All danumed® Gravity Feeding Sets are  
 • Latex-free • DEHP-free • BPA-free

## danumed® Gravity Universal

With **ENPlus Spike** for connection to danumed® FeedBags or pre-filled feeding bags and with additional **Universal Adapter** with vent valve for connection to crown caps and wide-neck bottles.

REF	Product description	Pieces/sales unit
DA 1051 0110	danumed® Gravity Universal	30 pieces

The danumed® Gravity Universal is individually sterile packed. 30 pieces in a shipping box.



The danumed® Universal Adapter with vent valve is also available as a separate product:

REF	Product description	Pieces/sales unit
DA 1051 8003	danumed® Universal Adapter	30 pieces

## Related products to danumed® Gravity Spike/Universal:



danumed® FeedBag



danumed® Universal Adapter



danumed® GastroTube



### danumed Medizintechnik GmbH – Individual solutions are our standard.

We are an innovative, German medical technology company with the highest quality standards, many years of expertise in enteral nutrition and clear customer focus. We also meet individual customer requests with flexibility and reliability. Personal, competent, collaborative.

danumed Medizintechnik GmbH • Stadthof 18 • 93059 Regensburg • Germany  
 Tel + 49 (0) 941/ 890 588-80 • Fax + 49 (0) 941/ 890 588-81  
 info@danumed.com • [www.danumed.com](http://www.danumed.com)

